MEDICAL DEVICES

T

0,

0

CompoGuard

16.10.13

Prese 500m

C

FRESENIUS KABI

CompoGuard®

MIXING/WEIGHING DEVICE COMBINING SAFETY WITH EFFICIENCY IN WHOLE BLOOD DONATION





The NEW CompoGuard®

MORE donations MORE flexibility



Increased number of donations per fully charged battery

- New high capacity Lithium battery technology
- Enhanced capability of approximately 70 donations per battery charge for mobile donations
- External charging station available
- Battery easily exchangeable

Increased data management flexibility

- Donation data now able to be downloaded via dedicated Fresenius Kabi USB stick

MORE benefits

One further step towards paperless donation

 A new software feature enables the direct donor verification via name, surname and date of birth




Mobile donation system


- Robust transport case for mobile donation

Comfortable working station

- The Gooseneck allows the operation of the CompoGuard[®] without bending down
- The large in-line tray and the tube clamp enable comfortable loading of the blood bag
- Transport case is also mobile and comfortable table during donation

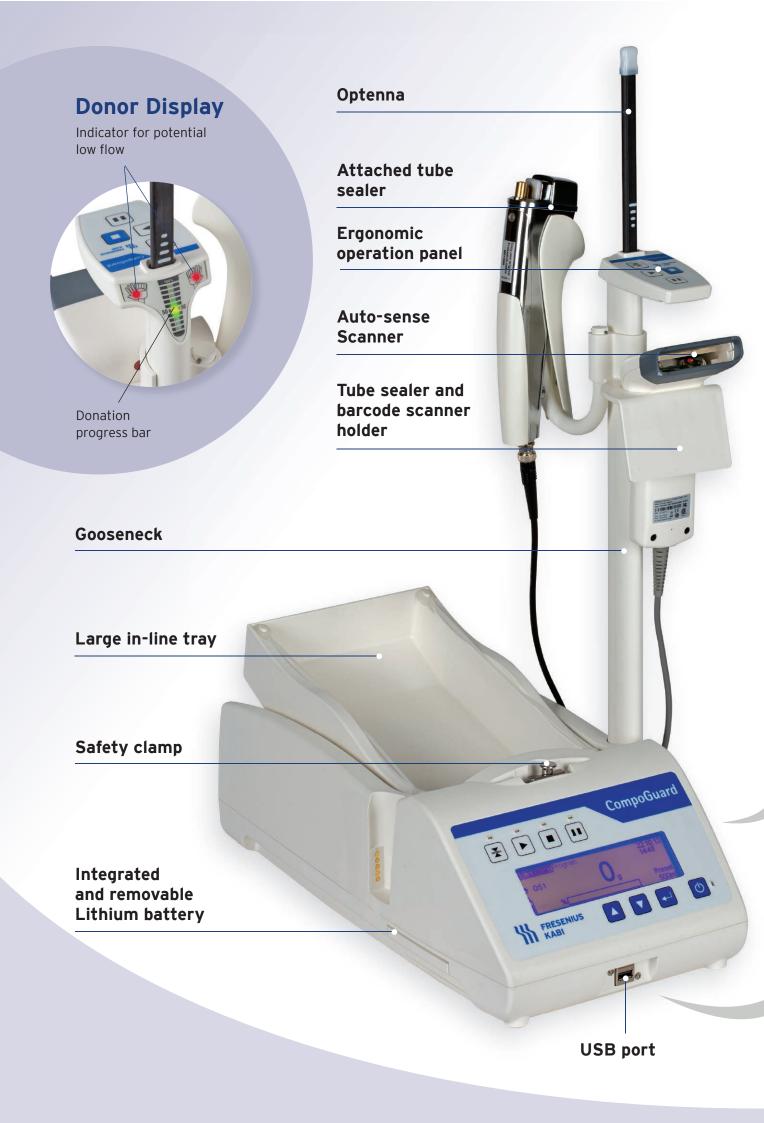

CompoMation Data Management System

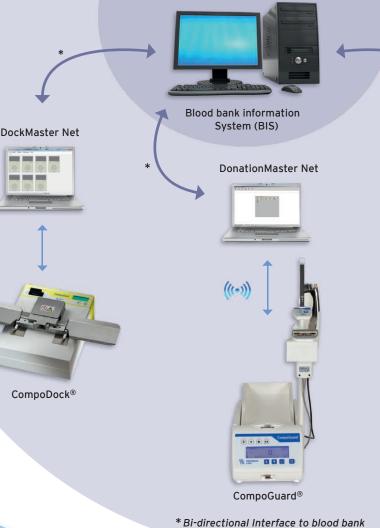

CompoMaster Net







CompoMat[®] G5



Fully

Data

Automated

- DonationMaster

Transmission

Net is part of the CompoMation Data Management System

Flexible data handling with DECT² technology

- Bi-directional data transfer
- Easy and flexible programming
- Reliable data recording, surveillance and documentation

² Digital Enhanced Cordless Telecommunications

Even more flexibility with USB assistance

- In addition to DECT and cable transmission, data transfer via USB stick is now possible
- Dedicated Fresenius Kabi USB stick

Trusted features of our CompoGuard[®] models



Optenna (optical alarm extension) for optimal donation supervision

- The donation event can be monitored from a distance
- The blinking Optenna leads you directly to the correct donor when action is needed



Safety clamp prevents air from entering the collection bag

- The donation line tube is fixed inside the safety clamp
- The tube clamp can only be opened after successful donation and sealing
- The heat sealing is integrated and documented by the CompoGuard®



Donor engagement to improve donation process

 Donors can directly monitor their donation progress



In-line tray for optimal results

- Integrated filter holder to support proper mixing
- Large tray for optimal bag positioning and accurate weighing

Blood Flow Indication The new CompoGuard®

can be programmed to indicate high and low flow events

CompoGuard[®] models



Ordering information

9029601	CompoGuard [®] Basic, 100-240 V
	Basic model

9029611 CompoGuard[®] Data, 100-240 V Basic model + DECT option + Gooseneck + USB port

9029621 CompoGuard[®] Complete, 100-240 V Basic model + DECT option + Gooseneck + USB port + sealing option + rechargeable battery

³ optional cable based data management

ORDERING INFORMATION

For more information such as literature, technical details and working procedures as well as for equipment please contact your local sales representative or us.

COMPOGUARD® MODELS

9029601	CompoGuard® Basic, 100 - 240 V Basic model (excl. battery)
9029611	CompoGuard® Data, 100 - 240 V
	Basic model + DECT option + Gooseneck
9029621	CompoGuard® Complete, 100 - 240 V
	Basic model + DECT option + sealing option + rechargeable battery + Gooseneck

DATA MANAGEMENT ACCESSORIES

M638421 M690041 9028891 9028341	Upgrade Option DECT for model "Basic" Upgrade Option DECT for model "Basic" (USA/Canada only) Barcode scanner with integrated auto-sense function DECT Fixedpart for radio transmission network, EU*
9029001	One per PC/network; includes power supply (230 V) and cable DECT Fixedpart clone, EU*
	Fixedpart clones all have the same DECT-ID number
9029361	DECT Fixedpart for radio transmission network, USA/CDN*
	One per PC/network; includes power supply (120 V) and cable
9029461	DECT Fixedpart clone, US A/CDN*
	Fixedpart clones all have the same DECT-ID number
M686791	USB adapter cable for Fixedpart
	With this adapter cable the Fixedpart can be connected to a PC with USB port
9029331	Donationmaster Net 2.5.x.x multilingual
9028391	Software licence for CompoGuard [®] Basic
	Key number for activating the data management; not necessary for models "Data" and "Complete"
9028331	Network converter RS485, EU (necessary for cable network only)*
	Includes power supply (230 V) and cable. One converter per PC/network
9028601	Network converter RS485, USA/CDN (necessary for cable network only)*
	Includes power supply (120 V) and cable. One converter per PC/network
M636701	Y-piece (necessary for cable network only)
	Number: one less than number of mixers in a network
M636671	Network cable, 5 m (necessary for cable network only)
	Number: one less than number of mixers in a network
9029541	USB 2.0 Stick 4GB for CompoGuard®

ACCESSORIES

9029471	Transport case for CompoGuard®
	Space for all accessories inside the transport case
9029481	Transport case with charging option
	External plug for in-case charging of battery
M684251	Power adapter with charging indication lamp 100 - 240V, 50 - 60 Hz*
M687051	Rechargeable battery for CompoGuard [®] , Li-Ion technology
9029491	4-fold charger for CompoGuard [®] , EU*
	For four batteries, 230 V
9029531	Gooseneck
	Includes holder for scanner and hand sealer; included in models "Data" and "Complete"
M657981	Disinfectable filter holder

* additional country/language versions available upon request

FRESENIUS KABI

caring for life

Fresenius Kabi AG

Else-Kröner-Straße 1 61352 Bad Homburg, Germany Phone: +49 (0) 61 72 686 5047 Fax: +49 (0) 61 72 686 4728 compoguard@fresenius-kabi.com www.fresenius-kabi.com The signs/names marked with \circledast are registered trademarks of the Fresenius Group in selected countries.

Your local sales organization: